

Description of iWhite Express set

iWhite Express is a revolutionary, dental safe teeth whitening device for immediate results. The ultra-thin sponge structure has a great absorption capacity and easily removes stubborn stains and dental plaque in hard-to-reach tooth surface areas. The sponge has to be activated with a few drops of the iWhite Express teeth whitening serum containing natural whitening enzymes, enamel strengthening actives and cooling agents for a breath freshening sensation. The combination of the sponge and serum leads to pearly white teeth and long term cleanliness of the tooth surface.



This ONE minute whitening treatment can be used once a week.

iWhite Express is an unique whitening treatment consisting of an ultra-thin sponge which needs to be activated by an innovative teeth whitening serum. iWhite Express is specifically designed to:


- **ACTIVELY REMOVE STUBBORN STAINS** caused by certain FOOD, DRINKS and TOBACCO from hard-to-reach tooth surface areas without harming the enamel. Thanks to the sponge technology, this innovative whitening system provides an immediate and superb whitening effect.
- **HELP PREVENTING** the formation of DENTAL PLAQUE and REMOVE STAINS that are bound to the dental plaque thanks to the presence of natural whitening enzymes in the iWhite Express teeth whitening serum.
- **STRENGTHEN** the enamel and **REMINERALIZE** the teeth by delivering highly available minerals such as Calcium & Magnesium.
- **NEUTRALIZES** mouth acids and **FIGHTS** sticky biofilm produced by bad bacteria in the mouth due to Xylitol.
- **GIVING you BREATH FRESHENING SENSATIONS** by using specialized cooling actives.

iWhite Express is clinically proven SAFE and EFFECTIVE.


Indications

iWhite Express was designed for dental cleaning and the removal of stubborn teeth stains. It also prevents plaque formation.

Product

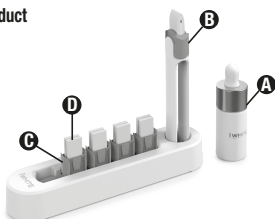

- A. iWhite Express teeth whitening serum
- B. iWhite Express teeth whitening holder
- C. iWhite Express plastic head
- D. iWhite Express stain eraser sponges

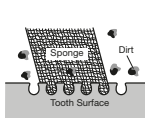


Ultra-thin, soft sponge based on nanotechnology with great absorption capacity to eliminate your stains and dental plaque.

Instructions for use

- Step 1:** Make sure that your hands are clean and dry
- Step 2:** Preparing for use:



- Step 3:** Gently shake the iWhite Express teeth whitening serum.

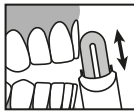

- Step 4:** Activate the stain eraser sponge with several drops of the iWhite Express teeth whitening serum and rub the teeth **GENTLY** with the flat side of the sponge and with vertical up and down movements. Make sure not to put too much pressure on the sponge. If the sponge gets damaged, use a new plastic head.
- Step 5:** Make sure the sponge stays moist during use, so if needed re-apply several drops of the iWhite Express teeth whitening serum.
- Step 6:** Rinse your mouth with water or a mouth-wash.

- Step 7:** Unclick the plastic head from the whitening holder. Remove the used stain eraser sponge from the plastic head and dispose.
- Step 8:** Wash the plastic head if you re-use it for the new stain eraser sponges from the refill bag.

How to use the iWhite Express stain eraser sponges from the refill bag?

- Step 1:** Make sure your hands are clean and dry.
- Step 2:** Place the plastic head with the curved side pointed downwards in the tray.
- Step 3:** Take a new sponge and fold it.
- Step 4:** Place the folded side of the sponge in the plastic head.



- Step 5:** Go further with step 2 of the instructions for use.

Maintenance

After each use, clean the holder with water. Do not use solvent-based cleaners.

Precautions

- Keep out of reach of children.
- Not recommended for children under 14 years.
- Not recommended for pregnant or breast-feeding women.
- Do not use more than directed.
- Do not perform this treatment in a prone position.
- This iWhite Express kit is for personal use only.
- Do not bite and/or chew the plastic parts and/or sponges.
- Avoid swallowing parts of the device, sponge and/or the iWhite Express teeth whitening serum.
- Rinse thoroughly with water in the event of contact of the iWhite Express teeth whitening serum with the eyes.
- Avoid direct contact of the plastic parts with the oral soft tissues (gums, tongue, lips) and the teeth.
- Damaged parts should no longer be used.
- Assure yourself the plastic head is correctly clicked on the holder.
- Always check if the sponge is correctly placed on the plastic heads, if not there is a risk the sponge comes loose and could be accidentally swallowed. If swallowed, contact your physician.
- Do not use the product when wearing braces, tooth retainers, veneers, tooth jewelry and/or piercings.
- Do not use this product if you suffer from gum problems. Discontinue the treatment if you experience pain or excessive sensitivity; if necessary contact your dental professional.
- Use of the product is inadvisable if you have allergic reactions to one of the ingredients of the iWhite Express teeth whitening serum and/or melamine.
- If you are undergoing treatment for any oral care condition, consult your dental professional prior to use.
- Do not use this product within 2 weeks in case of dental and/or orthodontic treatments (fillings, crowns, implants...).
- We recommend visiting your dentist every 6 months for healthy oral hygiene.
- Store in a cool dry place away from heat and direct sunlight.
- Keep away from dirt and dust.

Other information:

Manufacturer:
Sylphar nv
Xavier De Cocklaan 24
B-9831 Deurle, Belgium
Tel +32 9 321 70 80
www.iwhiteinstant.com

Content:

- 1 whitening holder
- 10 iWhite Express stain eraser sponges assembled on plastic heads
- 1 refill bag with 5 iWhite Express stain eraser sponges
- 1 iWhite Express teeth whitening serum (8 ml e)